



## **Manikin Integrated Data Acquisition (MIDAS) Initial Modifications**



By

**Nabih M. Alem**

**Aircrew Protection Division**

and

**James A. Lewis  
Robert M. Dillard**

**Research Support Division**

19950306 019

**December 1994**

Approved for public release; distribution unlimited.

**United States Army Aeromedical Research Laboratory  
Fort Rucker, Alabama 36362-0577**

## Notice

### Qualified requesters

Qualified requesters may obtain copies from the Defense Technical Information Center (DTIC), Cameron Station, Alexandria, Virginia 22314. Orders will be expedited if placed through the librarian or other person designated to request documents from DTIC.

### Change of address

Organizations receiving reports from the U.S. Army Aeromedical Research Laboratory on automatic mailing lists should confirm correct address when corresponding about laboratory reports.

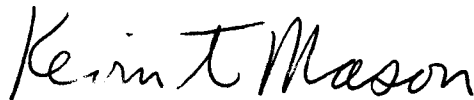
### Disposition

Destroy this document when it is no longer needed. Do not return it to the originator.

### Disclaimer

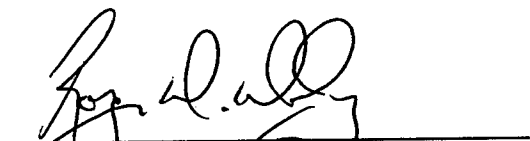
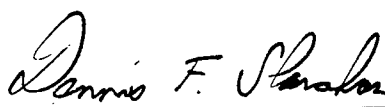
The views, opinions, and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy, or decision, unless so designated by other official documentation. Citation of trade names in this report does not constitute an official Department of the Army endorsement or approval of the use of such commercial items.

Reviewed:



KEVIN T. MASON  
LTC, MC, MFS  
Director, Aircrew Protection  
Division

Released for publication:

  
ROGER W. WILEY, O.D., Ph.D.  
Chairman, Scientific  
Review Committee  
DENNIS F. SHANAHAN  
Colonel, MC, MFS  
Commanding

REPORT DOCUMENTATION PAGE				Form Approved OMB No. 0704-0188	
1a. REPORT SECURITY CLASSIFICATION Unclassified			1b. RESTRICTIVE MARKINGS		
2a. SECURITY CLASSIFICATION AUTHORITY			3. DISTRIBUTION / AVAILABILITY OF REPORT Approved for public release, distribution unlimited		
2b. DECLASSIFICATION / DOWNGRADING SCHEDULE			5. MONITORING ORGANIZATION REPORT NUMBER(S)		
4. PERFORMING ORGANIZATION REPORT NUMBER(S) USAARL Report No. 95-11			7a. NAME OF MONITORING ORGANIZATION U.S. Army Medical Research and Materiel Command (Provisional)		
6a. NAME OF PERFORMING ORGANIZATION U.S. Army Aeromedical Research Laboratory		6b. OFFICE SYMBOL (If applicable) SGRD-UAD-IV	7b. ADDRESS (City, State, and ZIP Code) Fort Detrick Frederick, MD 21702-5012		
6c. ADDRESS (City, State, and ZIP Code) P.O. Box 620577 Fort Rucker, AL 36362-0577			9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER		
8a. NAME OF FUNDING / SPONSORING ORGANIZATION		8b. OFFICE SYMBOL (If applicable)	10. SOURCE OF FUNDING NUMBERS		
8c. ADDRESS (City, State, and ZIP Code)			PROGRAM ELEMENT NO. 62787A 30	PROJECT NO. 62787A878	TASK NO. ED
			WORK UNIT ACCESSION NO. 141		
11. TITLE (Include Security Classification) Manikin integrated data acquisition system (MIDAS) initial modifications					
12. PERSONAL AUTHOR(S) Nabih M. Alem, James A. Lewis, and Robert M. Dillard					
13a. TYPE OF REPORT Final		13b. TIME COVERED FROM _____ TO _____		14. DATE OF REPORT (Year, Month, Day) 1994 December	
15. PAGE COUNT 25					
16. SUPPLEMENTARY NOTATION					
17. COSATI CODES			18. SUBJECT TERMS (Continue on reverse if necessary and identify by block number)		
FIELD	GROUP	SUB-GROUP	prototype, manikin, crash testing, data acquisition, anthropomorphic test device, injury assessment		
15	05				
05	09				
19. ABSTRACT (Continue on reverse if necessary and identify by block number)					
<p>A prototype manikin with internal data acquisition system (MIDAS) has been developed and received by the US Army Aeromedical Research Laboratory. In addition to a novel design of the spinal column and pelvis of the Hybrid III automotive manikin, the new MIDAS includes a built-in signal conditioning and acquisition electronics. This report documents the initial modifications to the manikin and includes a description of the external software (MIDAS 3.0) for control, communication, and posttest downloading and analysis of the data.</p>					
20. DISTRIBUTION / AVAILABILITY OF ABSTRACT <input checked="" type="checkbox"/> UNCLASSIFIED/UNLIMITED <input type="checkbox"/> SAME AS RPT. <input type="checkbox"/> DTIC USERS			21. ABSTRACT SECURITY CLASSIFICATION Unclassified		
22a. NAME OF RESPONSIBLE INDIVIDUAL Chief, Science Support Center			22b. TELEPHONE (Include Area Code) 205-255-6907		22c. OFFICE SYMBOL SGRD-UAX-SI

## Contents

	Page
List of figures .....	2
Introduction .....	5
Materials and methods .....	6
Electronics hardware .....	6
Data acquisition settings .....	6
External software .....	7
Results .....	7
External software (MIDAS 3.0) .....	8
Discussion .....	9
Summary and conclusions .....	10
References .....	11
Appendix A. Manikin operating procedures .....	12
Appendix B. External manikin software, MIDAS version 3.0 .....	14

Accession For	
NTIS CRA&I	<input checked="" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification .....	
By .....	
Distribution /	
Availability Codes	
Dist	Avail and/or Special
A-1	

## List of figures

Figure		Page
B-1	Screen printout of the test identification menu which allows the user to enter a test label and additional comment associated with the test. ....	15
B-2.	Screen printout of the data extraction menu where the user selects the source and destination of the data to be extracted .....	15
B-3.	Screen printout of the data processing and analysis menu. This menu opens several submenus where various types of signal analyses are allowed .....	16
B-4.	Screen printout of the menu to setup the data acquisition sampling rate and the parameters for impact detection and data extraction from the credit card .....	16
B-5.	Screen printout of the menu to define channel labels and designate those active channels to be extracted after the test .....	17
B-6.	Screen printout of the menu to define the amplifiers gains in the data acquisition system, and to label the engineering units of each active channel. ....	17
B-7.	Screen printout of the menu to enter the transducers sensitivities of the manikin, as supplied by the manufacturers or determined from calibrations .....	18
B-8.	Screen printout of the menu to backup credit card data and other labels and settings for later retrieval and analysis .....	18
B-9.	Screen printout of the manikin controls menu which allows communication with the manikin internal data acquisition sub-system over serial communication port .....	19
B-10.	Screen printout of the manikin controls menu when the "monitor" mode is activated allowing synchronous serial communication with the manikin .....	19
B-11.	Screen printout of the menu to setup the colors and communication port parameters. Choosing colors becomes important for laptop LCD screens .....	20

# List of figures (Continued).

Figure		Page
B-12.	Screen printout of data processing submenu which allows the selection of up to six channels to be graphically displayed, in raw form, on the screen . . . . .	20
B-13.	Screen printout of the data processing submenu allowing the selection of up to four channels to be plotted in portrait format "A" hard copy for quick look. . . . .	21
B-14.	Screen printout of the data processing submenu allowing the selection of up to nine channels to be plotted in portrait format "B" hard copy for quick look . . . . .	21
B-15.	Screen printout of the data processing sub-menu allowing the selection of up to four channels to be plotted in landscape format "C" hard copy for quick look . . . . .	22
B-16.	Screen printout of the data processing submenu allowing the selection of up to nine channels to be plotted in landscape format "D" hard copy for quick look . . . . .	22
B-17.	Screen printout of the data processing submenu allowing the selection of up to eight channels to be plotted in two columns in landscape format "E" hard copy for quick look . . . . .	23
B-18.	Screen printout of the data processing submenu allowing the selection of up to 18 channels to be plotted in landscape format "F" hard copy for quick look . . . . .	23
B-19.	Screen printout of the data processing submenu showing test information and allowing the printing of hard copy of pertinent information . . . . .	24
B-20.	Screen printout of the data processing submenu showing (bottom third) parameters for integral/derivative analysis and allowing the user to modify the parameters . . . . .	24
B-21.	Screen printout of the data processing submenu showing parameters for over-plot of two signals. This analysis also allows plotting of a single channel . . . . .	25
B-22.	Screen printout of the data processing submenu showing (bottom third) the parameters selected to analyze a triaxial cluster of transducers . . . . .	25

---

---

This page left blank intentionally.

---

---

## Introduction

A prototype manikin with internal data acquisition system (MIDAS) has been developed and received by the U.S. Army Aeromedical Research Laboratory (USAARL), Fort Rucker, Alabama, for initial evaluations. New features in this anthropomorphic manikin include a novel design of the spinal column that offers more flexibility in compression and twisting than other manikins, a pelvis that contains the signal conditioning and data acquisition electronics required in most impact tests, and two new upper legs where rechargeable nickel-cadmium cells are placed to provide power supply for the instrumentation and data acquisition system.

As with most prototype systems, MIDAS must be evaluated before it can be put to full use in crash tests and other impact exposures. The evaluation process is simple: expose the manikin to impact tests, collect its biodynamic response signals, then compare the impact response to that of other validated manikins such as the original Hybrid III. This requires that the MIDAS internal electronics and data acquisition system operate properly and allow the data to be extracted in a short time.

Initial bench testing of MIDAS's signal conditioning electronics pointed out areas of potential problems which must be addressed prior to any sled or field testing. First, the battery packs which were installed in the manikin femurs were not supplying adequate power to the electronics. Additionally, the power supply that runs off the standard 120 VAC line and which was delivered with the prototype introduced unacceptable levels of noise to the signal amplifiers. Therefore, a new power supply was needed to replace the current battery packs and provide power from an external source through an umbilical line.

The second area where difficulties were encountered was the software for data acquisition, downloading, and display, which was delivered with the system. Since our immediate goal was the evaluation of the biodynamic response, it was essential that we are able to modify the internal settings of the electronics and to immediately examine the data captured after each test for quick assessment of the manikin performance. Our bench tests also pointed out the need for reliable and flexible software which allows us to achieve these two simple goals: modify the settings and examine the data as quickly as possible.

This report is intended to document the initial improvements made to the MIDAS prototype since its delivery. Because it was only a prototype, modifications were made only to bring the manikin to an acceptable level of reliability during sled tests.

## Materials and methods

### Electronics hardware

The external power supply was designed to convert 120 VAC to provide several regulated DC voltages to MIDAS's internal computer CPU, memory, disks, transducers, and amplifiers. The box delivered with the manikin had to be reworked by first removing the battery charger subsystem from the power supply box itself. Second, the soldering joints of the cables had to be redone because of potential breakage during anticipated sled mounting of the power supply box.

A serious problem that was identified after extensive debugging was the so-called "common mode problem." This problem caused the signal amplifiers to saturate due to the transducer excitation being tied to system ground. The problem was corrected by floating (relative to system ground) the 15-Volt excitation voltage of each transducer power supply.

Another shortcoming was cumbersome access to the analog signals from the 24 transducers installed in the manikin. Access to these signals is necessary for debugging and maintenance, for bypassing the data acquisition system in case of failure, and for recording the analog signals on an external data recording medium. At one stage during the development and testing of MIDAS, the contractor has built into the wiring access to all the transducer signals. Since this was not part of the contract and because of space constraints, the access cables were removed. In order for USAARL to make its evaluation, it was necessary to reconstruct the access cable and connect it to a breakout box where BNC-type connectors may be used. This proved later to be a valuable modification when it became necessary to record the signals external to the MIDAS hardware.

### Data acquisition settings

Precise measurements require accurate knowledge and control of the sampling rate and the gains for each amplifier in the manikin internal data acquisition system. Therefore, it was essential to verify the values of sampling rate and gains claimed by the contractor. We used the voltage insertion method, where a low level sine wave signal, generated by a precision calibrator, was used to replace a transducer and provide an input to the amplifier. By measuring the output of the amplifier, we determined its actual gain as the ratio of its output over the input, and its zero offset as the midpoint of the peak-to-peak range of the output.

Although it was possible to determine the amplifiers' gains by monitoring analog input and output signals, it was necessary to determine the sampling rate of the A/D convertor from the digitized data after it was stored on the thin credit-card memory, then dumped onto a PC disk file and examined. This was done after the first version of the MIDAS software became operational.

## External software

As delivered, the operation of MIDAS requires the user to turn on the power supply, modify the amplifier gains and other data acquisition parameters to accommodate specific test conditions, conduct the tests, then remove the credit-card memory module (where the manikin internal computer stored the data) for down-loading and analysis of the digitized data. A new MIDAS program was written to externally control the amplifier settings and sampling rate, and to extract the data from the credit-card and display it for quick look and assessment of the quality of the data acquisition.

Because of its intended use, the MIDAS program was designed to improve the communication by automating most of the functions which were required to be repeated by the user. The program was written in Fortran language (Microsoft Corp, 1988) to be run on an IBM-PC computer under a MS-DOS environment. All test data would be downloaded and archived along with descriptive and calibration data. These archived files may be accessed by MIDAS for further analysis and plotting at a later time after completing a test series. The program produces hard copies of raw and processed signals and summaries of analyzed signals. The program offered significant improvements in the operation of the manikin and, although several additions were contemplated, the first version was deemed adequate for evaluating the results from sled tests.

## Results

The bench evaluation of the electronics, wiring, and data acquisition system produced the following observations. Refer to Frisch, Boulay, and Alem, 1994, for a description of the delivered MIDAS and explanation of the components mentioned here.

1. The signal conditioner circuitry has no offset capability. Potentiometers on the board that are connected to the sample-and-hold (S/H) amplifiers have no effect on the direct analog or S/H output.
2. The voltage substitution calibration, as currently provided, is not useful. The resistor calibration (RCAL) is not implemented.
3. Gain and sample rate settings cannot be verified and sometimes are in error. Occasionally the gain sets to gain code 0, and at other times the sample rate sets to lowest value, as evidenced by a very dim data card light-emitting diode (LED).
4. For proper operation of the sensors, the power supplies required rewiring to "float" the grounds. The following supplies were affected: +15 volts for all of the accelerometers, +15 volts for all of the load cells, and +10 volts for the angular accelerometer in the head.
5. The specified rates for sampling are only approximate values.

6. Sampling at the highest two rates, 5 kilo-Hertz (KHz) and 10 KHz, is not possible with the current system which writes 48 channels to the data card.
7. The manikin's internal wiring harness is difficult to check or correct. It should have been of a more replaceable design, so that sensors or cables could be replaced easily. The copper tape shielding is ineffective.
8. The definitions of channels 18 and 20 appear to be swapped from their original designations.
9. Several signal conditioner hybrid circuits are bad, in one or both channels of the hybrid.
10. The L5 load cell channels have noise problems; so do the Head-X and the T1-Y channels.
11. Channel 19 is nonfunctional, producing a 147 Hertz noise signal.
12. Signal conditioner module card #1 has bad U-11 chip.

#### External software (MIDAS 3.0)

The MIDAS software initially was written to replace the software delivered as part of the contract because that software did not allow sufficient flexibility and reliability in operating the manikin, nor adequate speed in downloading and examining the acquired data. The latest version of MIDAS (version 3.0) was the result of actual tests and experience in the manikin's operation as intensive bench testing was conducted.

The step-by-step procedure to operate the manikin is summarized in Appendix A. This requires the MIDAS software and a PC connected through an RS-232 cable to the manikin internal computer. The capabilities of the MIDAS (version 3.0) software are summarized in Appendix B in the form of a series of screen printouts of the program menus. Although shown in black and white, these "screens" are produced on a VGA monitor in colors which are used for emphasis of key words or important text, and for demarcation between blocks of information. The following features are the highlights of the software.

1. The program consists of over 200 subroutines, written in Microsoft Fortran 5 and developed in its entirety at USAARL.
2. The program is interactive when necessary and automated when possible, and is entirely menu driven. All of the major submenus in MIDAS are reproduced in Appendix A as screen printouts.

3. The program requires an IBM-PC compatible running MS-DOS 5.0 or higher.
4. The PC hardware should have a VGA color monitor, although a laptop LCD screen may be used after adjusting the "color" settings to suit the particular screen.
5. No mouse is needed but a serial port (COM1 or COM2) must be available to communicate with the manikin internal computer.
6. MIDAS software is designed to run independent of the manikin to analyze, plot, and produce hard copies of the analyzed data after all the data from a series of tests have been archived.
7. Hard copies are produced in PCL/5 and HP-GL/2 languages and require printers capable of interpreting these languages, such as the HP Laserjet 3 or 4.

### Discussion

As with most prototype systems, some adjustments to the delivered version of the system was necessary. The initial modifications were designed primarily to put the manikin in an acceptable and reliable working condition so its performance as a test device could be evaluated. The reliability of the power supply, the calibration of the amplifiers, and the identification and repair (when possible) of known problems were of the utmost concern. Just as critical was the down-load and quick-look capability without which the manikin simply could not have been tested efficiently.

The manikin was designed to be self contained, i.e., capable of operation from an internal power source. The design of the prototype included two sets of nickel-cadmium (nicad) batteries installed in the femur in two cylindrical arrays. Although this approach was discussed and approved during the development phase of the manikin, the execution of this design by a subcontractor resulted in a product that functioned below our expectations. A remedy to this problem was set aside while more important issues were addressed. Thus, the manikin was to operate from external power supply until the internal battery supply issue has been resolved. Indeed, an external power supply was delivered but it also had problems: noise, ground loops, poor shielding, weak soldering connections and unnecessary wiring. All these problem areas were addressed and resolved in the more reliable configuration which was fabricated. This allowed us to advance to the next stage in this research program: the actual sled and controlled drop tests.

Equally important was our ability to program the internal amplifiers' gains and adjust the sampling rate of the manikin data acquisition system to meet the different impact test requirements. The manikin had default settings for all amplifiers and sampling rate which were stored on an erasable/programmable ROM chip in the MIDAS computer and read each time the manikin was turned on. Access to these stored defaults required the "undressing" of the manikin

and the removal of the CPU board, then its installation in a separate but equivalent computer for re-coding. The alternative to this time-consuming task was the transmittal of the desired gain values over the serial communication port every time the manikin was turned on to override the default values.

This approach was an acceptable solution if all goes well. Unfortunately, it took nearly half an hour of interaction between the user and the manikin internal software before the new settings had been transmitted. Frequently, the transmittal was so erratic because of poor COM port hardware and software drivers or because of unexplainable malfunction in the MIDAS internal electronics that it was necessary to reset the internal electronics by turning off the manikin power supply then switching it back on to restart the entire operation. This frustrating and inefficient process was one of the two primary reasons for writing our own external MIDAS control software.

The other motive for rewriting the MIDAS software was our inability to examine the data immediately after it was acquired. This essential requirement was somewhat met by the software as it was delivered. In fact, we were able to quickly remove the credit card and "download" the raw data to a PC file using utilities supplied by the manufacturer of the credit card. Unfortunately, plotting one impact pulse, which may be of 20 or 30 millisecond duration, required the manual plotting of the entire channel signal (1200 msec) in small segments and, when all segments were finally plotted, did not give the user sufficient timing and magnitude information to judge the quality of the test. The tedious process was lengthy and often was abandoned after examining one or two channels inconclusively.

Therefore, the MIDAS program was designed to truly assist the user in making informed decisions regarding a test. Just as the manikin could not have been operated without a reliable power supply, it would have been unwise to operate this new prototype without the benefit of "quick look" and "hard copy" capabilities.

### Summary and conclusions

As with most prototypes, initial modifications were necessary to the hardware and the software delivered with the manikin in order to bring it to an acceptable operating status. These include improvements to the external power supply to be used temporarily in lieu of the internal battery system, and a PC-based program (MIDAS) which facilitates the settings of internal amplifiers gains and sampling rate, and allows quick look at the down-loaded data immediately following a test. Eventually, it will be necessary to redesign the internal battery-based power supply so that field tests can be conducted independent of external connections. Also, it will be necessary to modify the internal data acquisition software to improve serial communication and to allow easier modification and retention of default data acquisition parameters.

### References

Frisch, P., Boulay, W., and Alem, N. 1994. Design and development of an enhanced biodynamic manikin. Phase I report for contract DAMD17-90-C-0116. Fort Rucker, AL: U.S. Army Aeromedical Research Laboratory. USAARL contractor report No. CR 94-1.

Microsoft Corporation. 1989. Fortran reference manual, version 5.0 for MS-DOS operating systems.

## Appendix A.

### MANIKIN operating procedures.

1. Connect RS-232 cable between manikin and PC communications port 1 or 2, as set in the USAARL manikin integrated data acquisition system (MIDAS) program.
2. Install the memory module (credit-card) into the manikin chassis slot.
3. Power-up manikin. The green LED on control unit indicates +5 volt power is present. Wait approximately 20 seconds - a short flash of the red LED indicates end of boot sequence.
4. On PC, run the MIDAS program. Verify sample rate in the "Sampling setup" function, gains and correction factors in "Amplifiers" function, and sensitivities in "Transducers" function.
5. Using "Identify test" function, enter the new test ID and other desired information.
6. Use the "Manikin controls" function to send gains and sample rate (press F1). If manikin will not respond, send "escape" characters (press F8) to clear the communications port. If this doesn't help, check for correct communication settings under the "Colors and comm." section, set port at 1200 baud, 8 bits, no parity, and 1 stop bit. Verify cabling is connected to correct port. Manikin uses upper communications port connector on manikin input/output (I/O) card, which is located inside the rib cage.
7. "Arm" manikin when ready (press F9). Observe that red LED is on and bright. A dim LED indicates improper sample rate. If LED goes off after 6 seconds (at 2500 Hz sample rate) or after 16 seconds (at 1000 Hz sample rate), then a false trigger has occurred. Turn off manikin power, ensure trigger signal is at ground potential (0 volts), and reapply power to manikin. Resend sample rate and gains after manikin boots up. A steady on LED indicates the proper "armed" condition (data is being written to the card in a circular buffer).
8. Disconnect the RS-232 cable from the manikin.
9. Run test. A 0-to-5 volt transition on the trigger signal input will end the 200-points/channel circular buffer and start the post-trigger data recording. When the red LED goes out, the 12,000 post-trigger data points/channel have been written to the card (storage complete).

10. Remove data card. Insert data card in the PC, perform "Backup card data" operation, saving the card image to hard disk. Use "Extract signals" to retrieve data from the data card (or from a saved file). Use the "Process" function to view signals.

11. Turn off manikin power if done, or press RESET and go to step 2 for the next test.

## Appendix B

External manikin software.

MIDAS version 3.0.

This appendix contains screen printouts of most menus and sub-menus which are available in the MIDAS program. Because they are printed in black and white, much of the color highlights used to emphasize various key words and text, and to separate blocks of information on the screen are, of course, missing. However, these "screens" do provide a printed documentation in lieu of actually running the MIDAS program itself.

Test: LX65-56		MIDAS-3.0	Manikin Integrated DAS	13:01:03
Identify test				
Extract signals	Data ID/Label: LX65-56			
Process signals	Description of data:			
Sampling setup	Dynamic testing of MIDAS prototype crash manikin			
Define channels	Sled acceleration pulse: 10.3 peak G, 23.1 mph			
Amplifiers	Seat pitch: +35 degrees			
Transducers				
Backup card data				
Manikin controls				
Colors and comm.				
Quit to DOS	Identify test data with a unique label to be used in naming binary and exchange data files. Descriptive text may also be added to identify test conditions.			

Figure B-1. Screen printout of the test identification menu which allows the user to enter a test label and additional comment associated with the test.

Test: LX65-56		MIDAS-3.0	Manikin Integrated DAS	13:06:32
Identify test	F1	EXTRACT data from THINCARD into RAM, and store them in MIDAS Internal Format on file	20 signals	LX65-56.MIF
▶ Extract signals	F2	SELECT from listed files ⇒ the name and type of MIDAS file to retrieve from disk	LX65-56.MIF LX65-55.MIF LX65-52.MIF LX65-45.MIF	
Process signals	F3	RETRIEVE the selected .MIF disk file into RAM memory.	LX65-44.MIF LX65-43.MIF	
Sampling setup	F4	RESTORE data from 3 backup files now on directory C:\NBDL\BACK\	LX65-56.ini LX65-56.pre LX65-56.dat	
Define channels	F5	EXPORT data to text file:	LX65-56.MXP	
Amplifiers				
Transducers				
Backup card data				
Manikin controls				
Colors and comm.				
Quit to DOS	Extract data from CC memory and save on hard disk in the desired format. Retrieve binary files from hard disk for plotting or export in another format.			

Figure B-2. Screen printout of the data extraction menu where the user selects the source and destination of the data to be extracted.

Test: LX65-56	MIDAS-3.0	Manikin Integrated DAS	13:06:57
Identify test	F1	Raw channels	
Extract signals	F2	Test/channel information	
► Process signals	F3	[3] Integral/derivative	
Sampling setup	F4	[2] Over-plot signals	
Define channels	F5	[3] Integral/derivative	
Amplifiers	F6	[6] Head triax & HIC	
Transducers	F7	[5] Triax and resultant	
Backup card data	F8	[5] Triax and resultant	
Manikin controls	F10	VGA: 640x480, 80x30	
Colors and comm.			
Quit to DOS			
Plot signals for screen quick look or hardcopy print. To print, LaserJet III printer port must already be defined with the DOS PRINT command.			

Figure B-3. Screen printout of the data processing and analysis menu. This menu opens several submenus where various types of signal analyses are allowed.

Test: LX65-56	MIDAS-3.0	Manikin Integrated DAS	13:07:29
Identify test	Sampling rate:	2500.0 Hz/channel (nominal)	
Extract signals	Actual:	2512.5 Hz per channel	
Process signals		.398 msec interval	
► Sampling setup	Sampling rate must match rate sent to manikin, other parameters may be re-defined after test.		
Define channels	Test duration:	597.0 msec = 1500 pts/ch	
Amplifiers	Trigger method:	SOFTWARE a signal FALLS BELOW threshold	
Transducers	Start of data:	-99.5 msec 250 pts BEFORE trigger	
Backup card data	Trigger signal:	Ch 15 = SLED	
Manikin controls	Threshold:	41 mV = 10.1 G BELOW initial level	
Colors and comm.			
Quit to DOS			
Define sampling rate, durations of pre- and post-test calibration lengths, pre-trigger delay and duration of test, all depending on the size of CC memory.			

Figure B-4. Screen printout of the menu to setup the data acquisition sampling rate and the parameters for impact detection and data extraction from the credit card.

Test: LX65-56	MIDAS-3.0	Manikin Integrated DAS	13:07:56
Identify test	(20)		A/D
Extract signals	ON	Label	channel
Process signals	01	AX-HEAD	Head X (forward) accel
Sampling setup	02	AY-HEAD	Head Y (lateral) accel
Define channels	03	AZ-HEAD	Head Z (vertical) accel
Amplifiers	04	AAZ-HEAD	Head Y (pitch) accel
Transducers	05	FZ-NECK	Neck Z (vertical) force
Backup card data	06	FX-NECK	Neck X (forward) force
Manikin controls	07	MY-NECK	Neck Y (pitch) moment
Colors and comm.	08	AX-T1	T1 X (forward) accel
Quit to DOS	09	AY-T1	T1 Y (lateral) accel
	10	AZ-T1	T1 Z (vertical) accel
	11	FZ-T1	T1 Z (vertical) force
	12	FX-T1	T1 X (forward) force
	13	MY-T1	T1 Y (pitch) moment
	14	MX-T1	T1 X (roll) moment
	15	SLED	Sled (external) accel
		Define and describe channels. Name engineering units for the signals. Designate signals which are to be extracted, stored, and displayed.	

Figure B-5. Screen printout of the menu to define channel labels and designate those active channels to be extracted after the test.

Test: LX65-56	MIDAS-3.0	Manikin Integrated DAS	13:09:21			
Identify test	(20)	DASS	Correction	Units @		
Extract signals	Ch	Label	Gain	Exp/actual	no load	Units
Process signals	01	AX-HEAD	100	1.000	.0000	G
Sampling setup	02	AY-HEAD	25	1.087	.0000	G
Define channels	03	AZ-HEAD	100	1.000	.0000	G
Amplifiers	04	AAZ-HEAD	100	1.000	.0000	krad/s2
Transducers	05	FZ-NECK	500	1.000	.0000	kN
Backup card data	06	FX-NECK	250	1.000	.0000	kN
Manikin controls	07	MY-NECK	250	1.000	.0000	kN-m
Colors and comm.	08	AX-T1	25	1.000	.0000	G
Quit to DOS	09	AY-T1	20	1.000	.0000	G
	10	AZ-T1	25	1.000	.0000	G
	11	FZ-T1	1000	1.000	.0000	kN
	12	FX-T1	250	1.000	.0000	kN
	13	MY-T1	250	1.000	.0000	kN-m
	14	MX-T1	250	1.000	.0000	kN-m
	15	SLED	1	.1000	.0000	G
Amplifier gains are defined here before sending them to the DASS in the manikin. Also, define slopes and intercepts of channel Out/In for correct calibration.						

Figure B-6. Screen printout of the menu to define the amplifiers gains in the data acquisition system, and to label the engineering units of each active channel.

Test: LX65-56	MIDAS-3.0	Manikin Integrated DAS	13:09:40		
Identify test	(20)	mV / Unit	Units/ADC		
Extract signals	Ch	Label	(@ gain=1)	Units	
Process signals	01	AX-HEAD	.5050	.2417E-01	G
	02	AY-HEAD	.5380	.9076E-01	G
	03	AZ-HEAD	.5570	.2192E-01	G
Sampling setup	04	AAY-HEAD	3.960	.3083E-02	krad/s2
	05	FZ-NECK	1.300	.1878E-02	kN
Define channels	06	FX-NECK	2.800	.1744E-02	kN
	07	MY-NECK	85.98	.5679E-04	kN-m
Amplifiers	08	AX-T1	2.430	.2009E-01	G
	09	AY-T1	2.440	.2501E-01	G
▶ Transducers	10	AZ-T1	2.440	.2001E-01	G
	11	FZ-T1	.8913	.1370E-02	kN
Backup card data	12	FX-T1	2.511	.1945E-02	kN
	13	MY-T1	47.60	.1026E-03	kN-m
Manikin controls	14	MX-T1	46.67	.1046E-03	kN-m
	15	SLED	4.070	.2999	G
Colors and comm.					
Quit to DOS	Transducers sensitivities (e.g., mV/G) may be edited. Overall factors to convert A/D counts to mechanical units are shown here for verification.				

Figure B-7. Screen printout of the menu to enter the transducers sensitivities of the manikin, as supplied by the manufacturers or determined from calibrations.

Test: LX65-56	MIDAS-3.0	Manikin Integrated DAS	13:10:43
Identify test	Page	Save credit card data on backup files	
Extract signals	Backup subdirectory:	C:\NBDL\BACK\	
Process signals	Configuration file:	1.INI	
Sampling setup	Pre-trigger data:	2.PRE	
Define channels	Post-trigger data:	3.DAT	
Amplifiers	Esc	Do NOT backup, CANCEL edit changes	
Transducers	F1	Do NOT backup, but SAVE file names	
Backup card data	F9	RE-WRITE modified file: LX65-56.MIF	
Manikin controls	Backup card data using TCREAD.EXE utility to create pre- & post-trigger (.PRE, .DAT) data files. Manikin configuration will also be saved on (.CFG) file.		
Colors and comm.			
Quit to DOS			

Figure B-8. Screen printout of the menu to backup credit card data and other labels and settings for later retrieval and analysis.

Test: LX65-56	MIDAS-3.0	Manikin Integrated DAS	13:11:03
Identify test	F1	send SAMPLING rate and ALL GAINS to DASS	
Extract signals	F2	send only SAMPLING rate = 2512 Hz to DASS	
Process signals	F3	send GAIN = 100x <7> to Ch 01 = AX-HEAD	
Sampling setup	F4	request pre-test calibration in DASS	
Define channels	F5	MONITOR: terminal mode over COM2 serial port	
Amplifiers	F6	request post-test calibration in DASS	
Transducers	F7	future use	
Backup card data	F8	send 3 ESC to DASS (to clear hung COM2 port)	
► Manikin controls	F9	arm manikin (put DASS in acquisition mode)	
Colors and comm.			
Quit to DOS		Control Panel for manikin: initializes manikin to specified configuration, arms manikin for trigger detection, and downloads acquired data from card.	

Figure B-9. Screen printout of the manikin controls menu which allows communication with the manikin internal data acquisition subsystem over serial communication port.

Test: LX65-56	MIDAS-3.0	Manikin Integrated DAS	13:12:08
Identify test	F1	send SAMPLING rate and ALL GAINS to DASS	
Extract signals	F2	send only SAMPLING rate = 2512 Hz to DASS	
Process signals	F3	send GAIN = 100x <7> to Ch 01 = AX-HEAD	
Sampling setup	F4	request pre-test calibration in DASS	
Define channels	F5	MONITOR: terminal mode over COM2 serial port	
Amplifiers			
Transducers			
Backup card data			
► Manikin controls			
Colors and comm.			
Quit to DOS			
			Port status byte COM2: 0110 0000
			Del      End clear    exit window   monitor
		Control Panel for manikin: initializes manikin to specified configuration, arms manikin for trigger detection, and downloads acquired data from card.	

Figure B-10. Screen printout of the manikin controls menu when the "monitor" mode is activated allowing synchronous serial communication with the manikin.

Test: LX65-56		MIDAS-3.0		Manikin Integrated DAS		13:12:32	
Identify test Extract signals Process signals Sampling setup Define channels Amplifiers Transducers Backup card data Manikin controls ▶ Colors and comm. Quit to DOS	Static	Lo	Hi	<div style="border: 1px solid black; padding: 5px; display: inline-block;"> <div style="text-align: center;">COLORS</div> <div>Static      Static</div> <div>Dynamic    Dynamic</div> <div>Key        Key      Key</div> <div>Message    Message</div> </div>			
	Dynamic	Lo	Hi				
	Buttons	Lo	Hi				
	Messages	Lo	Hi				
	F	B	colors				
Speaker		Port	Baud	Parity	Word	Stop	
On		COM2	1200	None	8-bit	1-bit	
Change color scheme of this program, and parameters of serial port that communicates with the manikin.							

Figure B-11. Screen printout of the menu to setup the colors and communication port parameters. Choosing colors becomes important for laptop LCD screens.

F1	Raw channels	F2	Test/channel information
F3	[1] Cross-plot signals	F4	[2] Over-plot signals
F5	[3] Integral/derivative	F6	[4] 1st/2nd integrals
F7	[5] Triax and resultant	F8	[6] Head triax & HIC
F9	Execute all ◀▶ functions	F10	VGA: 640x480, 80x30
Fn or Enter screen output		Alt - Fn parameters	
		Tab form/graph	
		Esc exit	
01 = AX-HEAD	= Head X (forward) accel	01	<div style="border: 1px solid black; padding: 10px; display: inline-block;">           -----            -----            -----            -----            -----            -----            -----            -----            -----            -----         </div>
02 = AY-HEAD	= Head Y (lateral) accel	02	
03 = AZ-HEAD	= Head Z (vertical) accel	03	
04 = AAY-HEAD	= Head Y (pitch) accel	05	
05 = FZ-NECK	= Neck Z (vertical) force	06	
06 = FX-NECK	= Neck X (forward) force	15	
07 = MY-NECK	= Neck Y (pitch) moment		
08 = AX-T1	= T1 X (forward) accel		
09 = AY-T1	= T1 Y (lateral) accel		
10 = AZ-T1	= T1 Z (vertical) accel		
		Analysis span <span style="border-bottom: 1px solid black; display: inline-block; width: 100px;"></span> 6APRTST3	

Figure B-12. Screen printout of data processing sub-menu which allows the selection of up to six channels to be graphically displayed, in raw form, on the screen.

F1	Raw channels	F2	Test/channel information
F3	[1] Cross-plot signals	F4	[2] Over-plot signals
F5	[3] Integral/derivative	F6	[4] 1st/2nd integrals
F7	[5] Triax and resultant	F8	[6] Head triax & HIC
F9	Execute all ◀▶ functions	F10	VGA: 640x480, 80x30
Ctrl - Fn hardcopy output		Alt - Fn parameters	Tab form/graph
			Esc exit
01 = AX-HEAD	= Head X (forward) accel	01	-----
02 = AY-HEAD	= Head Y (lateral) accel	02	-----
03 = AZ-HEAD	= Head Z (vertical) accel	03	-----
04 = AAY-HEAD	= Head Y (pitch) accel	04	-----
05 = FZ-NECK	= Neck Z (vertical) force		
06 = FX-NECK	= Neck X (forward) force		
07 = MY-NECK	= Neck Y (pitch) moment		
08 = AX-T1	= T1 X (forward) accel		
09 = AY-T1	= T1 Y (lateral) accel		
10 = AZ-T1	= T1 Z (vertical) accel		
		HARDCOPY	6APRTST3 ↑A

Figure B-13. Screen printout of the data processing submenu allowing the selection of up to four channels to be plotted in portrait format "A" hard copy for quick look.

F1	Raw channels	F2	Test/channel information
F3	[1] Cross-plot signals	F4	[2] Over-plot signals
F5	[3] Integral/derivative	F6	[4] 1st/2nd integrals
F7	[5] Triax and resultant	F8	[6] Head triax & HIC
F9	Execute all ◀▶ functions	F10	VGA: 640x480, 80x30
Ctrl - Fn hardcopy output		Alt - Fn parameters	Tab form/graph
			Esc exit
01 = AX-HEAD	= Head X (forward) accel	01	-----
02 = AY-HEAD	= Head Y (lateral) accel	02	-----
03 = AZ-HEAD	= Head Z (vertical) accel	03	-----
04 = AAY-HEAD	= Head Y (pitch) accel	04	-----
05 = FZ-NECK	= Neck Z (vertical) force	05	-----
06 = FX-NECK	= Neck X (forward) force	06	-----
07 = MY-NECK	= Neck Y (pitch) moment	07	-----
08 = AX-T1	= T1 X (forward) accel	08	-----
09 = AY-T1	= T1 Y (lateral) accel	09	-----
10 = AZ-T1	= T1 Z (vertical) accel		
		HARDCOPY	6APRTST3 2B

Figure B-14. Screen printout of the data processing submenu allowing the selection of up to nine channels to be plotted in portrait format "B" hard copy for quick look.

F1	Raw channels	F2	Test/channel information
F3	[1] Cross-plot signals	F4	[2] Over-plot signals
F5	[3] Integral/derivative	F6	[4] 1st/2nd integrals
F7	[5] Triax and resultant	F8	[6] Head triax & HIC
F9	Execute all ◀▶ functions	F10	VGA: 640x480, 80x30
Ctrl - Fn		Alt - Fn	parameters
hardcopy output		Tab	form/graph
		Esc	exit
01 = AX-HEAD	= Head X (forward) accel	01	-----
02 = AY-HEAD	= Head Y (lateral) accel	02	-----
03 = AZ-HEAD	= Head Z (vertical) accel	03	-----
04 = AAY-HEAD	= Head Y (pitch) accel	04	-----
05 = FZ-NECK	= Neck Z (vertical) force		
06 = FX-NECK	= Neck X (forward) force		
07 = MY-NECK	= Neck Y (pitch) moment		
08 = AX-T1	= T1 X (forward) accel		
09 = AY-T1	= T1 Y (lateral) accel		
10 = AZ-T1	= T1 Z (vertical) accel	HARDCOPY	6APRTST3 ↑C

Figure B-15. Screen printout of the data processing submenu allowing the selection of up to four channels to be plotted in landscape format "C" hard copy for quick look.

F1	Raw channels	F2	Test/channel information
F3	[1] Cross-plot signals	F4	[2] Over-plot signals
F5	[3] Integral/derivative	F6	[4] 1st/2nd integrals
F7	[5] Triax and resultant	F8	[6] Head triax & HIC
F9	Execute all ▶◀ functions	F10	VGA: 640x480, 80x30
Ctrl - Fn		Alt - Fn	parameters
hardcopy output		Tab	form/graph
		Esc	exit
01 = AX-HEAD	= Head X (forward) accel	01	-----
02 = AY-HEAD	= Head Y (lateral) accel	02	-----
03 = AZ-HEAD	= Head Z (vertical) accel	03	-----
04 = AAY-HEAD	= Head Y (pitch) accel	04	-----
05 = FZ-NECK	= Neck Z (vertical) force	05	-----
06 = FX-NECK	= Neck X (forward) force	06	-----
07 = MY-NECK	= Neck Y (pitch) moment	07	-----
08 = AX-T1	= T1 X (forward) accel	08	-----
09 = AY-T1	= T1 Y (lateral) accel	09	-----
10 = AZ-T1	= T1 Z (vertical) accel	HARDCOPY	6APRTST3 2D

Figure B-16. Screen printout of the data processing submenu allowing the selection of up to nine channels to be plotted in landscape format "D" hard copy for quick look.

F1	Raw channels	F2	Test/channel information
F3	[1] Cross-plot signals	F4	[2] Over-plot signals
F5	[3] Integral/derivative	F6	[4] 1st/2nd integrals
F7	[5] Triax and resultant	F8	[6] Head triax & HIC
F9	Execute all ◀▶ functions	F10	VGA: 640x480, 80x30
Ctrl - Fn hardcopy output		Alt - Fn parameters	Tab form/graph
			Esc exit
01 = AX-HEAD	= Head X (forward) accel	01 -----	05 -----
02 = AY-HEAD	= Head Y (lateral) accel	02 -----	06 -----
03 = AZ-HEAD	= Head Z (vertical) accel	03 -----	07 -----
04 = AAY-HEAD	= Head Y (pitch) accel	04 -----	08 -----
05 = FZ-NECK	= Neck Z (vertical) force		
06 = FX-NECK	= Neck X (forward) force		
07 = MY-NECK	= Neck Y (pitch) moment		
08 = AX-T1	= T1 X (forward) accel		
09 = AY-T1	= T1 Y (lateral) accel		
10 = AZ-T1	= T1 Z (vertical) accel	HARDCOPY	6APRTST3 3E

Figure B-17. Screen printout of the data processing submenu allowing the selection of up to eight channels to be plotted in two columns in landscape format "E" hard copy for quick look.

F1	Raw channels	F2	Test/channel information
F3	[1] Cross-plot signals	F4	[2] Over-plot signals
F5	[3] Integral/derivative	F6	[4] 1st/2nd integrals
F7	[5] Triax and resultant	F8	[6] Head triax & HIC
F9	Execute all ◀▶ functions	F10	VGA: 640x480, 80x30
Ctrl - Fn hardcopy output		Alt - Fn parameters	Tab form/graph
			Esc exit
01 = AX-HEAD	= Head X (forward) accel	01 -----	10 -----
02 = AY-HEAD	= Head Y (lateral) accel	02 -----	11 -----
03 = AZ-HEAD	= Head Z (vertical) accel	03 -----	12 -----
04 = AAY-HEAD	= Head Y (pitch) accel	04 -----	13 -----
05 = FZ-NECK	= Neck Z (vertical) force	05 -----	14 -----
06 = FX-NECK	= Neck X (forward) force	06 -----	15 -----
07 = MY-NECK	= Neck Y (pitch) moment	07 -----	16 -----
08 = AX-T1	= T1 X (forward) accel	08 -----	17 -----
09 = AY-T1	= T1 Y (lateral) accel	09 -----	18 -----
10 = AZ-T1	= T1 Z (vertical) accel	HARDCOPY	6APRTST3 1F

Figure B-18. Screen printout of the data processing sub-menu allowing the selection of up to 18 channels to be plotted in landscape format "F" hard copy for quick look.

F1	Raw channels	◄►	F2	Test/channel information
F3	[3] Integral/derivative		F4	[2] Over-plot signals
F5	[3] Integral/derivative		F6	[2] Over-plot signals
F7	[1] Cross-plot signals		F8	[2] Over-plot signals
F9	Execute all ◄► functions		F10	VGA: 640x480, 80x30
Ctrl - Fn hardcopy output      Ins mark for hardcopy      Esc exit				
Test LX65-56      1500 pts @ 2512.5 Hz = 597.0 ms long Descript.: Dynamic testing of MIDAS prototype crash manikin Sled acceleration pulse: 10.3 peak G, 23.1 mph Seat pitch: +35 degrees  20 channels: 01 - 20      Analysis: 59.7 ms Span: .0 to 59.7 ms.				

Figure B-19. Screen printout of the data processing submenu showing test information and allowing the printing of hard copy of pertinent information.

F1	Raw channels	F2	Test/channel information
F3	[3] Integral/derivative	F4	[2] Over-plot signals
◄► F5	[3] Integral/derivative	F6	[2] Over-plot signals
F7	[1] Cross-plot signals	F8	[2] Over-plot signals
F9	Execute all ◄► functions	F10	VGA: 640x480, 80x30
Fn or Enter screen output      Alt - Fn parameters      Tab form/graph			
Ctrl - Fn hardcopy output      Ins mark for hardcopy      Esc exit			
Signals: 05 FZ-NECK      First integral times 1.0000      First derivative times 1.0000 Filter: 100 Hz      Set integral BEG. = 0 (LP) 24 dB/oct Labels: Head axial (Z) force      Head Z (G) (kN) LX65-56 Caption: 7 Gx run, WITHOUT cushion			

Figure B-20. Screen printout of the data processing submenu showing (bottom third) parameters for integral/derivative analysis and allowing the user to modify the parameters.

F1	Raw channels	F2	Test/channel information
F3	[3] Integral/derivative	◄► F4	[2] Over-plot signals
F5	[3] Integral/derivative	F6	[2] Over-plot signals
F7	[1] Cross-plot signals	F8	[2] Over-plot signals
F9	Execute all ◄► functions	F10	VGA: 640x480, 80x30
Fn or Enter screen output      Alt - Fn parameters      Tab form/graph Ctrl - Fn hardcopy output      Ins mark for hardcopy      Esc exit			
Signals: 10 AZ-T1      (Left vertical)      00 Seat pit      (Right vertical) Filter: 100 Hz (LP) 24 dB/oct Labels: Neck (C1) axial force (kN) LX65-56 Caption: 7 Gx run WITHOUT cushion			

Figure B-21. Screen printout of the data processing submenu showing parameters for over-plot of two signals. This analysis also allows plotting of a single channel.

F1	Raw channels	F2	Test/channel information
F3	[3] Integral/derivative	F4	[2] Over-plot signals
F5	[3] Integral/derivative	F6	[2] Over-plot signals
◄► F7	[5] Triax and resultant	F8	[2] Over-plot signals
F9	Execute all ◄► functions	F10	VGA: 640x480, 80x30
Fn or Enter screen output      Alt - Fn parameters      Tab form/graph Ctrl - Fn hardcopy output      Ins mark for hardcopy      Esc exit			
Signals: 01 AX-HEAD      02 AY-HEAD      03 AZ-HEAD      Tri-axial resultant Filter: 100 Hz (LP) 24 dB/oct Labels: Head X      Head Y      Head Z      Head resultant accel (G)      (G)      (G)      (G) LX65-56 Caption: 7 Gx run WITHOUT cushion			

Figure B-22. Screen printout of the data processing submenu showing (bottom third) the parameters selected to analyze a triaxial cluster of transducers.

## Initial distribution

Commander, U.S. Army Natick Research,  
Development and Engineering Center  
ATTN: SATNC-MIL (Documents  
Librarian)  
Natick, MA 01760-5040

Chairman  
National Transportation Safety Board  
800 Independence Avenue, S.W.  
Washington, DC 20594

Commander  
10th Medical Laboratory  
ATTN: Audiologist  
APO New York 09180

Naval Air Development Center  
Technical Information Division  
Technical Support Detachment  
Warminster, PA 18974

Commanding Officer, Naval Medical  
Research and Development Command  
National Naval Medical Center  
Bethesda, MD 20814-5044

Deputy Director, Defense Research  
and Engineering  
ATTN: Military Assistant  
for Medical and Life Sciences  
Washington, DC 20301-3080

Commander, U.S. Army Research  
Institute of Environmental Medicine  
Natick, MA 01760

Library  
Naval Submarine Medical Research Lab  
Box 900, Naval Sub Base  
Groton, CT 06349-5900

Executive Director, U.S. Army Human  
Research and Engineering Directorate  
ATTN: Technical Library  
Aberdeen Proving Ground, MD 21005

Commander  
Man-Machine Integration System  
Code 602  
Naval Air Development Center  
Warminster, PA 18974

Commander  
Naval Air Development Center  
ATTN: Code 602-B  
Warminster, PA 18974

Commanding Officer  
Armstrong Laboratory  
Wright-Patterson  
Air Force Base, OH 45433-6573

Director  
Army Audiology and Speech Center  
Walter Reed Army Medical Center  
Washington, DC 20307-5001

Commander/Director  
U.S. Army Combat Surveillance  
and Target Acquisition Lab  
ATTN: SFAE-IEW-JS  
Fort Monmouth, NJ 07703-5305

Director  
Federal Aviation Administration  
FAA Technical Center  
Atlantic City, NJ 08405

Director  
Walter Reed Army Institute of Research  
Washington, DC 20307-5100

Commander, U.S. Army Test  
and Evaluation Command  
Directorate for Test and Evaluation  
ATTN: AMSTE-TA-M (Human Factors  
Group)  
Aberdeen Proving Ground,  
MD 21005-5055

Naval Air Systems Command  
Technical Air Library 950D  
Room 278, Jefferson Plaza II  
Department of the Navy  
Washington, DC 20361

Director  
U.S. Army Ballistic  
Research Laboratory  
ATTN: DRXBR-OD-ST Tech Reports  
Aberdeen Proving Ground, MD 21005

Commander  
U.S. Army Medical Research  
Institute of Chemical Defense  
ATTN: SGRD-UV-AO  
Aberdeen Proving Ground,  
MD 21010-5425

Commander  
USAMRMC  
ATTN: SGRD-RMS  
Fort Detrick, Frederick, MD 21702-5012

HQ DA (DASG-PSP-O)  
5109 Leesburg Pike  
Falls Church, VA 22041-3258

Harry Diamond Laboratories  
ATTN: Technical Information Branch  
2800 Powder Mill Road  
Adelphi, MD 20783-1197

U.S. Army Materiel Systems  
Analysis Agency  
ATTN: AMXSU-PA (Reports Processing)  
Aberdeen Proving Ground  
MD 21005-5071

U.S. Army Ordnance Center  
and School Library  
Simpson Hall, Building 3071  
Aberdeen Proving Ground, MD 21005

U.S. Army Environmental  
Hygiene Agency  
ATTN: HSHB-MO-A  
Aberdeen Proving Ground, MD 21010

Technical Library Chemical Research  
and Development Center  
Aberdeen Proving Ground, MD  
21010-5423

Commander  
U.S. Army Medical Research  
Institute of Infectious Disease  
ATTN: SGRD-UIZ-C  
Fort Detrick, Frederick, MD 21702

Director, Biological  
Sciences Division  
Office of Naval Research  
600 North Quincy Street  
Arlington, VA 22217

Commandant  
U.S. Army Aviation  
Logistics School ATTN: ATSQ-TDN  
Fort Eustis, VA 23604

Headquarters (ATMD)  
U.S. Army Training  
and Doctrine Command  
ATTN: ATBO-M  
Fort Monroe, VA 23651

IAF Liaison Officer for Safety  
USAF Safety Agency/SEFF  
9750 Avenue G, SE  
Kirtland Air Force Base  
NM 87117-5671

Naval Aerospace Medical  
Institute Library  
Building 1953, Code 03L  
Pensacola, FL 32508-5600

Command Surgeon  
HQ USCENTCOM (CCSG)  
U.S. Central Command  
MacDill Air Force Base, FL 33608

Director  
Directorate of Combat Developments  
ATTN: ATZQ-CD  
Building 515  
Fort Rucker, AL 36362

U.S. Air Force Institute  
of Technology (AFIT/LDEE)  
Building 640, Area B  
Wright-Patterson  
Air Force Base, OH 45433

Henry L. Taylor  
Director, Institute of Aviation  
University of Illinois-Willard Airport  
Savoy, IL 61874

Chief, National Guard Bureau  
ATTN: NGB-ARS  
Arlington Hall Station  
111 South George Mason Drive  
Arlington, VA 22204-1382

AAMRL/HEX  
Wright-Patterson  
Air Force Base, OH 45433

Commander  
U.S. Army Aviation and Troop Command  
ATTN: AMSAT-R-ES  
4300 Goodfellow Bouvelard  
St. Louis, MO 63120-1798

U.S. Army Aviation and Troop Command  
Library and Information Center Branch  
ATTN: AMSAV-DIL  
4300 Goodfellow Boulevard  
St. Louis, MO 63120

Federal Aviation Administration  
Civil Aeromedical Institute  
Library AAM-400A  
P.O. Box 25082  
Oklahoma City, OK 73125

Commander  
U.S. Army Medical Department  
and School  
ATTN: Library  
Fort Sam Houston, TX 78234

Commander  
U.S. Army Institute of Surgical Research  
ATTN: SGRD-USM  
Fort Sam Houston, TX 78234-6200

Air University Library  
(AUL/LSE)  
Maxwell Air Force Base, AL 36112

Product Manager  
Aviation Life Support Equipment  
ATTN: SFAE-AV-LSE  
4300 Goodfellow Boulevard  
St. Louis, MO 63120-1798

Commander and Director  
USAE Waterways Experiment Station  
ATTN: CEWES-IM-MI-R,  
CD Department  
3909 Halls Ferry Road  
Vicksburg, MS 39180-6199

Commanding Officer  
Naval Biodynamics Laboratory  
P.O. Box 24907  
New Orleans, LA 70189-0407

Assistant Commandant  
U.S. Army Field Artillery School  
ATTN: Morris Swott Technical Library  
Fort Sill, OK 73503-0312

Mr. Peter Seib  
Human Engineering Crew Station  
Box 266  
Westland Helicopters Limited  
Yeovil, Somerset BA20 2YB UK

U.S. Army Dugway Proving Ground  
Technical Library, Building 5330  
Dugway, UT 84022

U.S. Army Yuma Proving Ground  
Technical Library  
Yuma, AZ 85364

AFFTC Technical Library  
6510 TW/TSTL  
Edwards Air Force Base,  
CA 93523-5000

Commander  
Code 3431  
Naval Weapons Center  
China Lake, CA 93555

Aeromechanics Laboratory  
U.S. Army Research and Technical Labs  
Ames Research Center, M/S 215-1  
Moffett Field, CA 94035

Sixth U.S. Army  
ATTN: SMA  
Presidio of San Francisco, CA 94129

Commander  
U.S. Army Aeromedical Center  
Fort Rucker, AL 36362

Strughold Aeromedical Library  
Document Service Section  
2511 Kennedy Circle  
Brooks Air Force Base, TX 78235-5122

Dr. Diane Damos  
Department of Human Factors  
ISSM, USC  
Los Angeles, CA 90089-0021

U.S. Army White Sands  
Missile Range  
ATTN: STEWS-IM-ST  
White Sands Missile Range, NM 88002

Director, Airworthiness Qualification Test  
Directorate (ATTC)  
ATTN: STEAT-AQ-O-TR (Tech Lib)  
75 North Flightline Road  
Edwards Air Force Base, CA 93523-6100

Ms. Sandra G. Hart  
Ames Research Center  
MS 262-3  
Moffett Field, CA 94035

Commander  
USAMRMC  
ATTN: SGRD-UMZ  
Fort Detrick, Frederick, MD 21702-5009

Commander  
U.S. Army Health Services Command  
ATTN: HSOP-SO  
Fort Sam Houston, TX 78234-6000

U. S. Army Research Institute  
Aviation R&D Activity  
ATTN: PERI-IR  
Fort Rucker, AL 36362

Commander  
U.S. Army Safety Center  
Fort Rucker, AL 36362

U.S. Army Aircraft Development  
Test Activity  
ATTN: STEBG-MP-P  
Cairns Army Air Field  
Fort Rucker, AL 36362

Commander  
USAMRMC  
ATTN: SGRD-PLC (COL R. Gifford)  
Fort Detrick, Frederick, MD 21702

TRADOC Aviation LO  
Unit 21551, Box A-209-A  
APO AE 09777

Netherlands Army Liaison Office  
Building 602  
Fort Rucker, AL 36362

British Army Liaison Office  
Building 602  
Fort Rucker, AL 36362

Italian Army Liaison Office  
Building 602  
Fort Rucker, AL 36362

Directorate of Training Development  
Building 502  
Fort Rucker, AL 36362

Chief  
USAHEL/USAAVNC Field Office  
P. O. Box 716  
Fort Rucker, AL 36362-5349

Commander, U.S. Army Aviation Center  
and Fort Rucker  
ATTN: ATZQ-CG  
Fort Rucker, AL 36362

Dr. Sehchang Hah  
Dept. of Behavior Sciences and  
Leadership, Building 601, Room 281  
U. S. Military Academy  
West Point, NY 10996-1784

Canadian Army Liaison Office  
Building 602  
Fort Rucker, AL 36362

German Army Liaison Office  
Building 602  
Fort Rucker, AL 36362

French Army Liaison Office  
USAAVNC (Building 602)  
Fort Rucker, AL 36362-5021

Australian Army Liaison Office  
Building 602  
Fort Rucker, AL 36362

Dr. Garrison Rapmund  
6 Burning Tree Court  
Bethesda, MD 20817

Commandant, Royal Air Force  
Institute of Aviation Medicine  
Farnborough, Hampshire GU14 6SZ UK

Defense Technical Information  
Cameron Station, Building 5  
Alexandria, VA 22304-6145

Commander, U.S. Army Foreign Science  
and Technology Center  
AIFRTA (Davis)  
220 7th Street, NE  
Charlottesville, VA 22901-5396

Commander  
Applied Technology Laboratory  
USARTL-ATCOM  
ATTN: Library, Building 401  
Fort Eustis, VA 23604

Commander, U.S. Air Force  
Development Test Center  
101 West D Avenue, Suite 117  
Eglin Air Force Base, FL 32542-5495

Aviation Medicine Clinic  
TMC #22, SAAF  
Fort Bragg, NC 28305

Dr. H. Dix Christensen  
Bio-Medical Science Building, Room 753  
Post Office Box 26901  
Oklahoma City, OK 73190

Commander, U.S. Army Missile  
Command  
Redstone Scientific Information Center  
ATTN: AMSMI-RD-CS-R  
/ILL Documents  
Redstone Arsenal, AL 35898

Aerospace Medicine Team  
HQ ACC/SGST3  
162 Dodd Boulevard, Suite 100  
Langley Air Force Base,  
VA 23665-1995

U.S. Army Research and Technology  
Laboratories (AVSCOM)  
Propulsion Laboratory MS 302-2  
NASA Lewis Research Center  
Cleveland, OH 44135

Commander  
USAMRMC  
ATTN: SGRD-ZC (COL John F. Glenn)  
Fort Detrick, Frederick, MD 21702-5012

Dr. Eugene S. Channing  
166 Baughman's Lane  
Frederick, MD 21702-4083

U.S. Army Medical Department  
and School  
USAMRDALC Liaison  
ATTN: HSMC-FR  
Fort Sam Houston, TX 78234

NVESD  
AMSEL-RD-NV-ASID-PST  
(Attn: Trang Bui)  
10221 Burbeck Road  
Fort Belvoir, VA 22060-5806

CA Av Med  
HQ DAAC  
Middle Wallop  
Stockbridge, Hants S020 8DY UK

Dr. Christine Schlichting  
Behavioral Sciences Department  
Box 900, NAVUBASE NLON  
Groton, CT 06349-5900

Commander  
Aviation Applied Technology Directorate  
ATTN: AMSAT-R-TV  
Fort Eustis, VA 23604-5577

COL Yehezkel G. Caine, MD  
Surgeon General, Israel Air Force  
Aeromedical Center Library  
P. O. Box 02166 I.D.F.  
Israel

HQ ACC/DOHP  
205 Dodd Boulevard, Suite 101  
Langley Air Force Base,  
VA 23665-2789

41st Rescue Squadron  
41st RQS/SG  
940 Range Road  
Patrick Air Force Base,  
FL 32925-5001

48th Rescue Squadron  
48th RQS/SG  
801 Dezonía Road  
Holloman Air Force Base,  
NM 88330-7715

HQ, AFOMA  
ATTN: SGPA (Aerospace Medicine)  
Bolling Air Force Base,  
Washington, DC 20332-6128

ARNG Readiness Center  
ATTN: NGB-AVN-OP  
Arlington Hall Station  
111 South George Mason Drive  
Arlington, VA 22204-1382

35th Fighter Wing  
35th FW/SG  
PSC 1013  
APO AE 09725-2055

66th Rescue Squadron  
66th RQS/SG  
4345 Tyndall Avenue  
Nellis Air Force Base, NV 89191-6076

71st Rescue Squadron  
71st RQS/SG  
1139 Redstone Road  
Patrick Air Force Base,  
FL 32925-5000

Director  
Aviation Research, Development  
and Engineering Center  
ATTN: AMSAT-R-Z  
4300 Goodfellow Boulevard  
St. Louis, MO 63120-1798

Commander  
USAMRMC  
ATTN: SGRD-ZB (COL C. Fred Tyner)  
Fort Detrick, Frederick, MD 21702-5012

Commandant  
U.S. Army Command and General Staff  
College  
ATTN: ATZL-SWS-L  
Fort Leavenworth, KS 66027-6900

ARNG Readiness Center  
ATTN: NGB-AVN-OP  
Arlington Hall Station  
111 South George Mason Drive  
Arlington, VA 22204-1382

Director  
Army Personnel Research Establishment  
Farnborough, Hants GU14 6SZ UK

Dr. A. Kornfield  
895 Head Street  
San Francisco, CA 94132-2813

ARNG Readiness Center  
AATN: NGB-AVN-OP  
Arlington Hall Station  
111 South George Mason Drive  
Arlington, VA 22204-1382

Cdr, PERSCOM  
ATTN: TAPC-PLA  
200 Stovall Street, Rm 3N25  
Alexandria, VA 22332-0413

HQ, AFOMA  
ATTN; SGPA (Aerospace Medicine)  
Bolling Air Force Base,  
Washington, DC 20332-6188